

Model No.:
BW121
Ssmart FamilyFit Scale

App Spec

1.1	Change model name and number	Kevin Yu	18/12/2013
1.0	Art work update	Kevin Yu	26-11-2013
0	First Release	Kevin Yu	18-09-2013
Rev.	Description	Issue	Date

Product Feature:


- One scale for different user group
- Measure weight and BMI
- BLE connection

App Features:


- Interface for input user data, e.g. age, height etc
- Separate different user group: Normal, baby, pregnancy, kid.
- Normal: Talking weight and extra large fonts for elderly
- Baby: growth weight chart
- Pregnancy: weight gain tracker
- Kid: height chart


System Requirement:

- iOS 6.0 or above
- Android 4.3 or above
- Support BLE

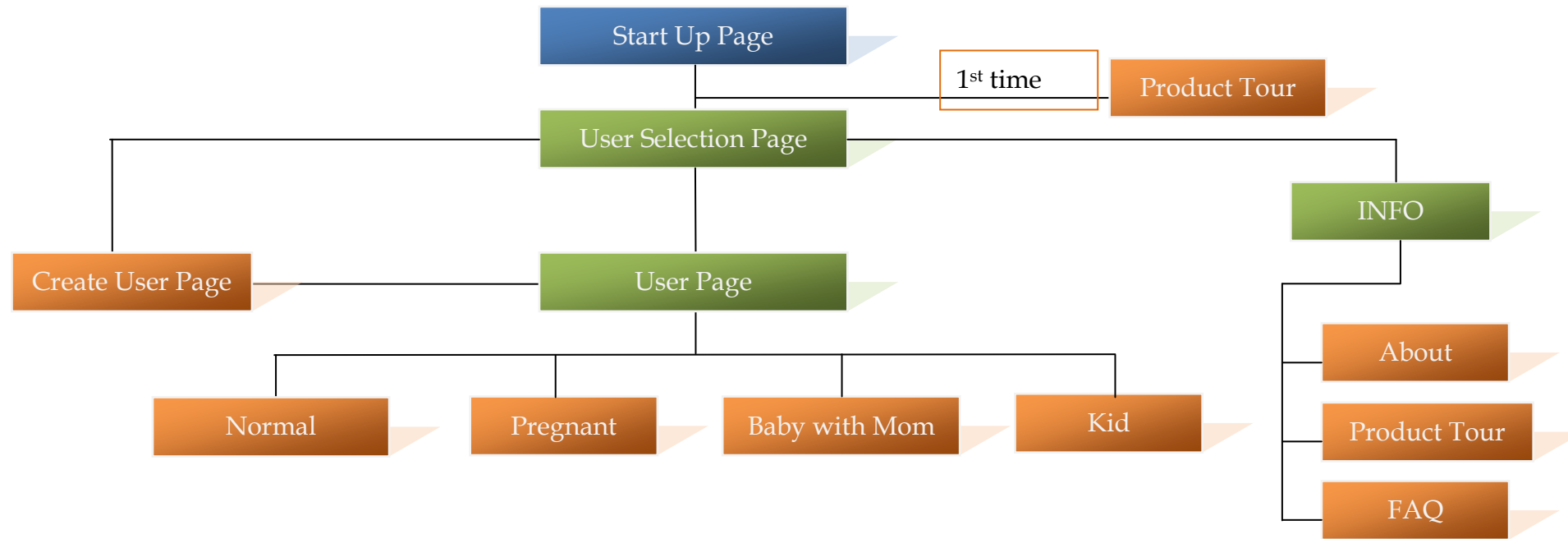

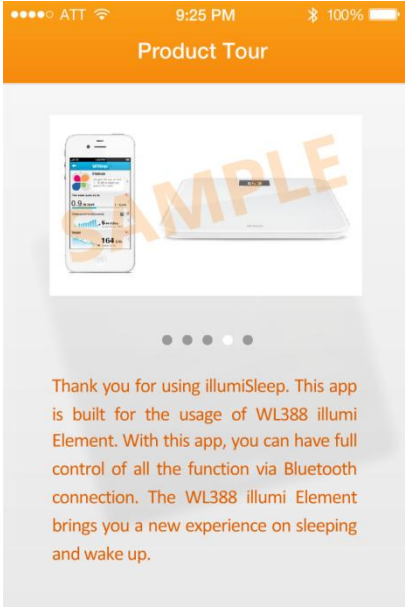



Screen Shot	Content	Remark
	<p>Startup Page</p> <ul style="list-style-type: none"> ➤ Show OS Logo and Product Shot ➤ Auto connect to main unit and show the connecting icon. TBC ➤ User still can move into User Selection Page although no BT connection. 	
	<p>Product Tour</p> <ul style="list-style-type: none"> ➤ The content of product tour will be defined in this file later 	<p>Reference WL388</p>

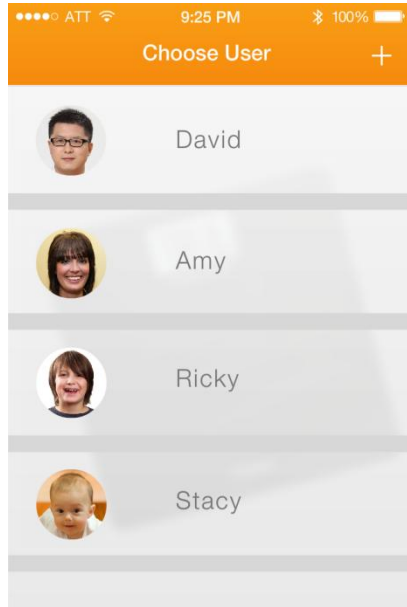

User Selection Page

- Show the thumbnail of all existing user a/c for selection
- Max number of user a/c is 10.
- A button for creating new user a/c
- An indication showing the connection status, i.e. connected, or disconnected

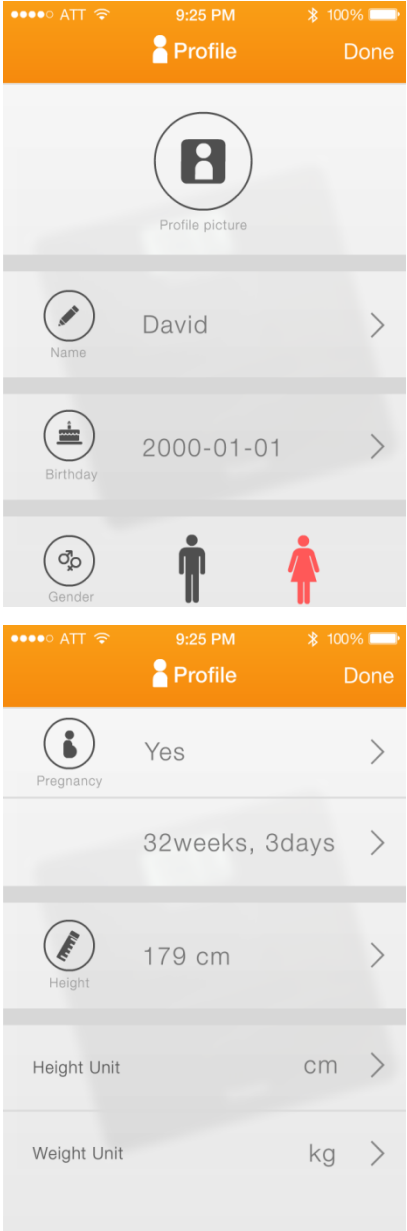Screen Shot	Content	Remark
 <p>The image contains two screenshots of a mobile application's 'Create User Page'. Both screenshots have an orange header bar with a 'Profile' icon and a 'Done' button. The status bar at the top shows 'ATT', signal strength, time '9:25 PM', and 100% battery.</p> <p>The top screenshot displays the following fields:</p> <ul style="list-style-type: none"> Profile picture: A circular icon with a person silhouette. Name: A text field containing 'David' with a pencil icon and a right arrow. Birthday: A date picker field showing '2000-01-01' with a calendar icon and a right arrow. Gender: A section with a gender selection icon (a circle with a cross) and two human icons, one male (black) and one female (red). <p>The bottom screenshot displays the following fields:</p> <ul style="list-style-type: none"> Pregnancy: A toggle switch field with a person icon and the text 'Yes' and a right arrow. Height: A text field showing '32weeks, 3days' with a ruler icon and a right arrow. Height Unit: A dropdown menu showing '179 cm' and 'cm' with a right arrow. Weight Unit: A dropdown menu showing 'kg' and a right arrow. 	<p>Create User Page</p> <ul style="list-style-type: none"> ➤ Input name ➤ Input photo ➤ Input date of birth ➤ Input gender <ul style="list-style-type: none"> ✧ If female, ask pregnancy and week ➤ Input height ➤ Select unit kg/lb, cm/ft 	

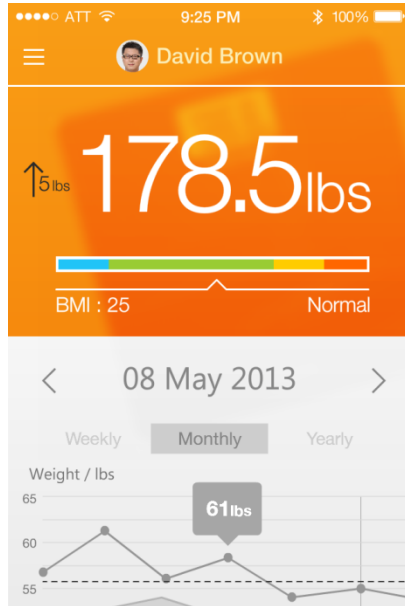

User Page

The User Page shows the measurement and history of user weighting info. Based on user info, they will be assigned to 4 different age groups, which are Normal, Pregnancy, Baby With Mom, Kid.


Normal (age 11-120)

- Show user name
 - Press on thumbnail or name to edit user info
- Show current weight (scale reading)
- Show current BMI
- Show weight chart
- Show BMI chart
- Show height chart


All chart can change by day/week/month/yr

$$\text{BMI} = \text{kg/m}^2$$

Classification	BMI(kg/m ²)	
	Principal cut-off points	Additional cut-off points
Underweight	<18.50	<18.50
Severe thinness	<16.00	<16.00
Moderate thinness	16.00 - 16.99	16.00 - 16.99
Mild thinness	17.00 - 18.49	17.00 - 18.49
Normal range	18.50 - 24.99	18.50 - 22.99
		23.00 - 24.99
Overweight	≥25.00	≥25.00
Pre-obese	25.00 - 29.99	25.00 - 27.49
		27.50 - 29.99


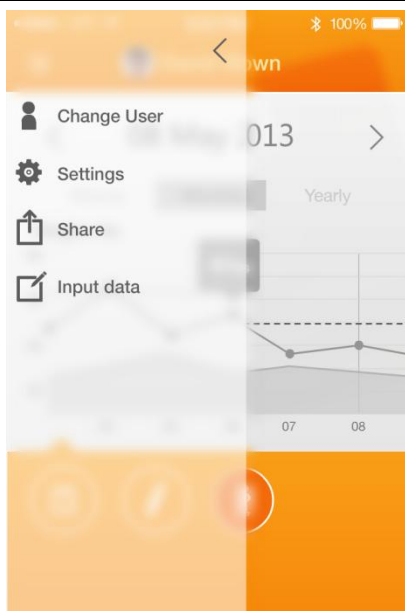


- Manual bar for option
 - Change user
 - Setting
 - Share to community
 - Manual input data

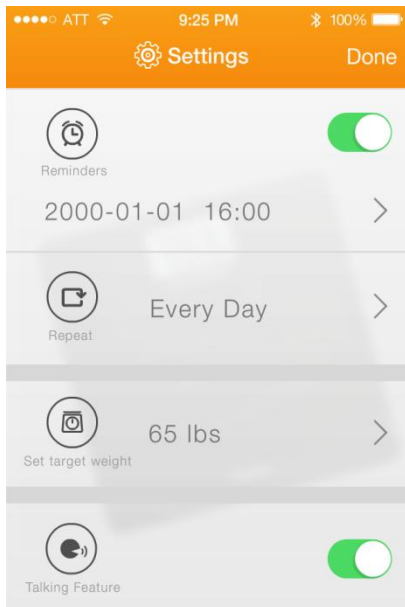


- Setting
 - ✧ Reminder on/off
 - ✧ Set alarm date and time for remind next measurement
 - ✧ Set target weight
 - ✧ Talking feature on/off
 - ✧ Auto Save
 - ✧ Product tour


Obese	≥ 30.00	≥ 30.00
Obese class I	30.00 - 34.99	30.00 - 32.49
		32.50 - 34.99
Obese class II	35.00 - 39.99	35.00 - 37.49
		37.50 - 39.99
Obese class III	≥ 40.00	≥ 40.00

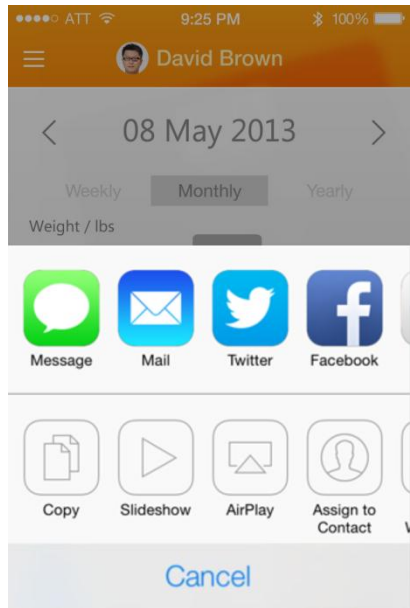
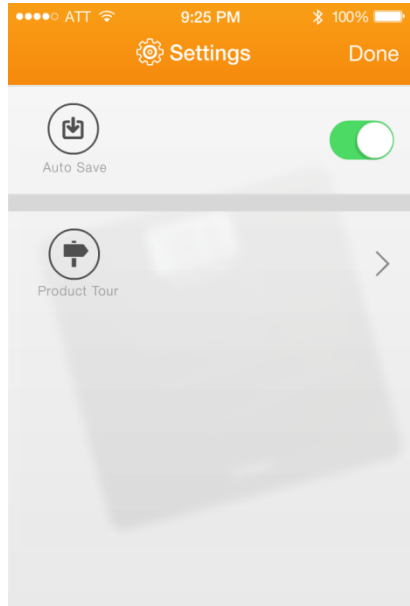


➤ Share to community

Only the weight chart can be shared

Screen Shot	Content	Remark
 <p>The screenshot displays a mobile application interface for pregnancy weight management. At the top, the status bar shows 'ATT' and '9:25 PM'. The app header includes a menu icon, a profile picture, and the name 'Amy'. The main display area has a red background and shows a large weight reading of '178.5 lbs' with an upward arrow and '5 lbs' indicating a change. Below the weight is a horizontal bar chart showing a BMI of 25 and a 'Normal' status. The bottom section features a date selector for '08 May 2013' and a line chart labeled 'BMI' with 'Weekly', 'Monthly', and 'Yearly' tabs. The 'Monthly' tab is selected, showing a line graph with a peak at 61.</p>	<p>Pregnancy (age 16-50)</p> <ul style="list-style-type: none"> ➤ Show user name <ul style="list-style-type: none"> - Press on thumbnail or name to edit user info ➤ Show current weight (scale reading) ➤ Show weight chart ➤ Show pregnancy weight gain reference (TBC) 	

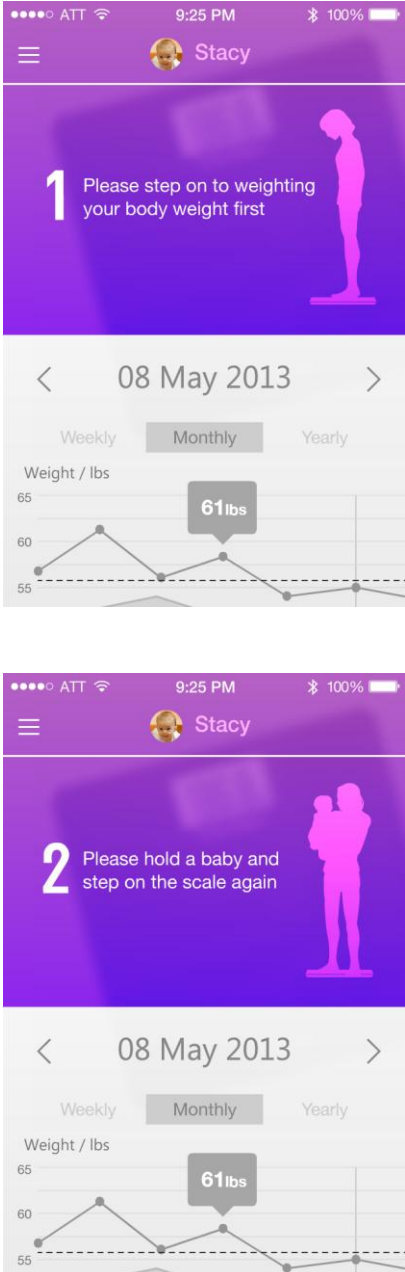Screen Shot	Content	Remark
	<p>Baby With Mom (age 0-1)</p> <ul style="list-style-type: none"> ➤ Step by step showing process for baby weight measurement ➤ Show baby name <ul style="list-style-type: none"> - Press on thumbnail or name to edit user info ➤ Show current weight (scale reading) ➤ Show weight chart, (with reference information, TBC) ➤ Show height chart, (with reference information, TBC) 	<p>Process for baby weight measurement:</p> <ul style="list-style-type: none"> ➤ Ask the mother (only) to stand on the scale ➤ Ask her to step down ➤ Ask the mother stand on the scale with her baby <p>All chart can change by day/week/month/yr</p> <p>Only the weight chart and height chart can be shared</p>

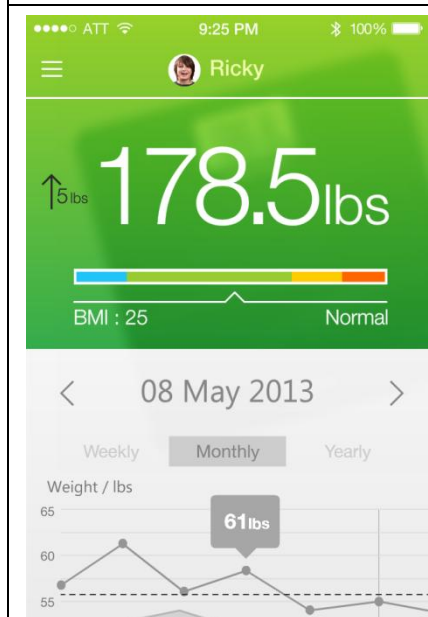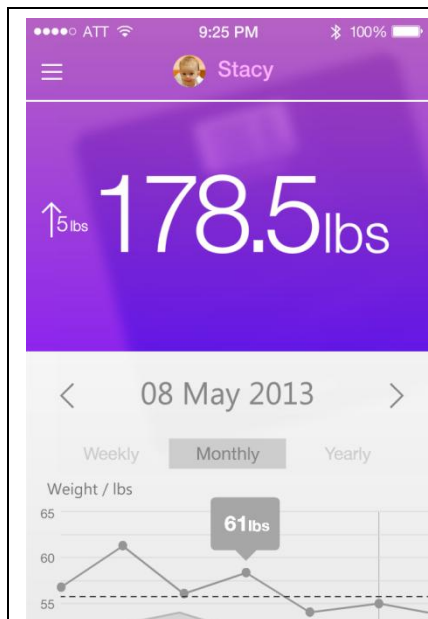

Kid (age 1-10)

- Show name
 - Press on thumbnail or name to edit user info
- Show current weight (scale reading)
- Show weight chart, (with reference information, TBC)
- Show height chart, (with reference information, TBC)
- Show BMI chart

All chart can change by day/week/month/yr

Only the weight chart, height chart and BMI chart can be shared

Screen Shot	Content	Remark
	Rotating the smart phone introduces a full screen reading and chart for better visualization	